

Monthly AED Check

External Indicators

Yes / No

- ☐ ☐ Status indicator is green
- ☐ ☐ Pads not expired. Battery not expired.
If expired replace. Order replacement battery
2 months before expiration date

Internal Indicators

Yes / No

- ☐ ☐ Adequate charge.
- ☐ ☐ Service indicator is off
- ☐ ☐ Pad connection indicator is off


Describe finding(s) if answered 'no' to any of the above questions. Reference user manual as needed for any abnormal indicators.

Additional Notes/Details

AED location: _____

Inspector name: _____

Inspection date: _____

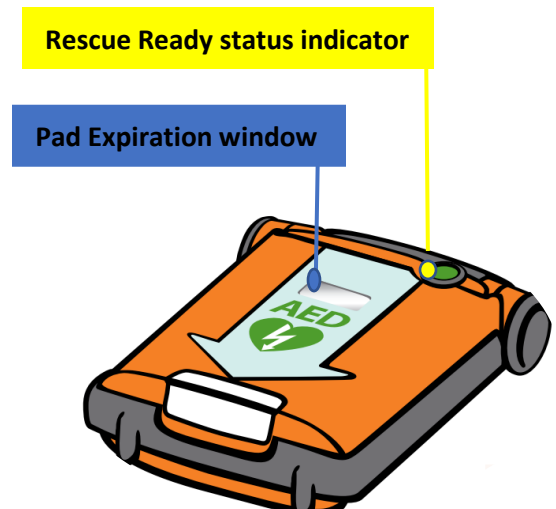

Battery capacity indicator

LEDs show the remaining battery capacity: 100%, 75%, 50%, 25%, 0% (red only).



Service Indicator indicates

AED requires service by authorized service personnel.



Defibrillation Pads indicator Indicates that the pads are incorrectly connected or unusable. Check the connection with the AED; check the placement and attachment to the patient. If connections are correct, replace pads.